

ardo melia single

WEARABLE
ELECTRIC BREASTPUMP



Convincing quality – fair price

Ardo Melia: a well-designed solution for occasional pumping – powerful, discreet and affordable.



The advantages

- **Mobile**
directly in the bra – without cables and tubes – for maximum of movement
- **Discreet**
worn in the bra, compact in design, quiet in operation
- **Relaxed**
hands-free, no need to constantly hold the pumpsets
- **Perfect fit**
 - Includes 2 breast shell/insert sizes: 19 and 24 mm
 - 13, 15, 17, 21 and 27 mm additionally available
 - Contains a strap extension for the nursing bra
- **Hospital-quality pumping**
Individually adjustable with four different modes, each with 15 vacuum steps – for maximum comfort and optimal milk flow
 - max. vacuum: 350 mbar/260 mmHG
 - max. cycles: 135/min



Mobile



Discreet



Relaxed

Pump hands-free – stay mobile

Proven technology – now also available as a single breast pump

The Ardo Melia Single was developed for mothers who occasionally express milk – wherever they are.

It helps to make expression comfortable, efficient and discreet. At the same time, the Ardo Melia Single offers sophisticated flexibility: it can be extended with a second pump if the need for pumping increases – so that mothers always have the right solution.

Individually adjustable – 4 modes & 15 vacuum steps

The Ardo Melia Single adapts perfectly to every mother: with four modes and 15 vacuum steps, pumping can be tailored to personal needs and individual preferences. These versatile settings ensure that every mother can quickly find her ideal combination of comfort and efficiency. The result: gentle, careful pumping that is also powerful and can be seamlessly integrated into everyday life – whether at home or on the go.



Comfortable fit

Thanks to the extensive range of accessories, the Melia Single can be easily adjusted – even for small nipple sizes. This adaptability not only ensures greater comfort, but also minimises discomfort. Mothers can therefore rely on a reliable, comfortable pumping experience.



Quiet operation – discreet in any situation

The pump operates quietly, allowing for discreet pumping in various environments, even in a shared space with other people.


Price-conscious & high quality

The Ardo Melia Single combines high-quality workmanship, functionality and practical accessories at a reasonable price. The Ardo Melia Single is a cost-effective choice as a practical solution for occasional pumping, without compromising on quality, comfort or efficiency.


Four different modes for maximum comfort and optimal milk flow

- Massage or stimulation mode gets the milk flowing


Massage: fast cycle, deep vacuum (vacuum manually adjustable)

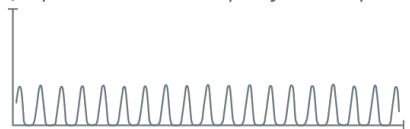

Stimulation: fast cycle, a vacuum peak every 10th time (vacuum manually adjustable)

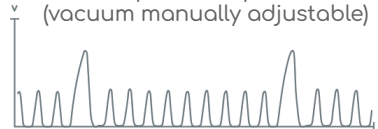


- As soon as the milk flows, switch to expression mode

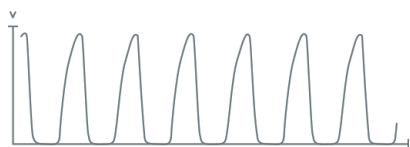

Slow cycle, high, comfortable vacuum (vacuum manually adjustable)

- Or alternatively automatic mode from the start.

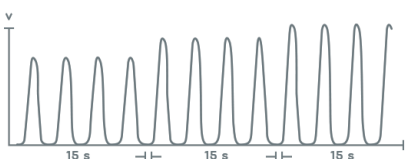

Starts at the lowest vacuum level and moves up one level every 15 seconds up to level 12 out of a total of 15 possible levels. Use the + or - button to switch to manual mode at any time and adjust the vacuum to your own requirements.

More facts



Powerful battery (1200 mAh) for short charging time (2 h) and long expression time (min. 100 min)



With strap extension



Silent



Soft breast shell for a snug fit and a comfortable feel



Free from BPA & BPS – for the baby's health



Closed system, everything stays hygienically clean



Large capacity of 180 ml / 6 fl. oz



Ardo is WHO code compliant




We from the Swiss family business Ardo medical guarantee the highest quality and first-class customer service.



instruction video



Simple and easy to use

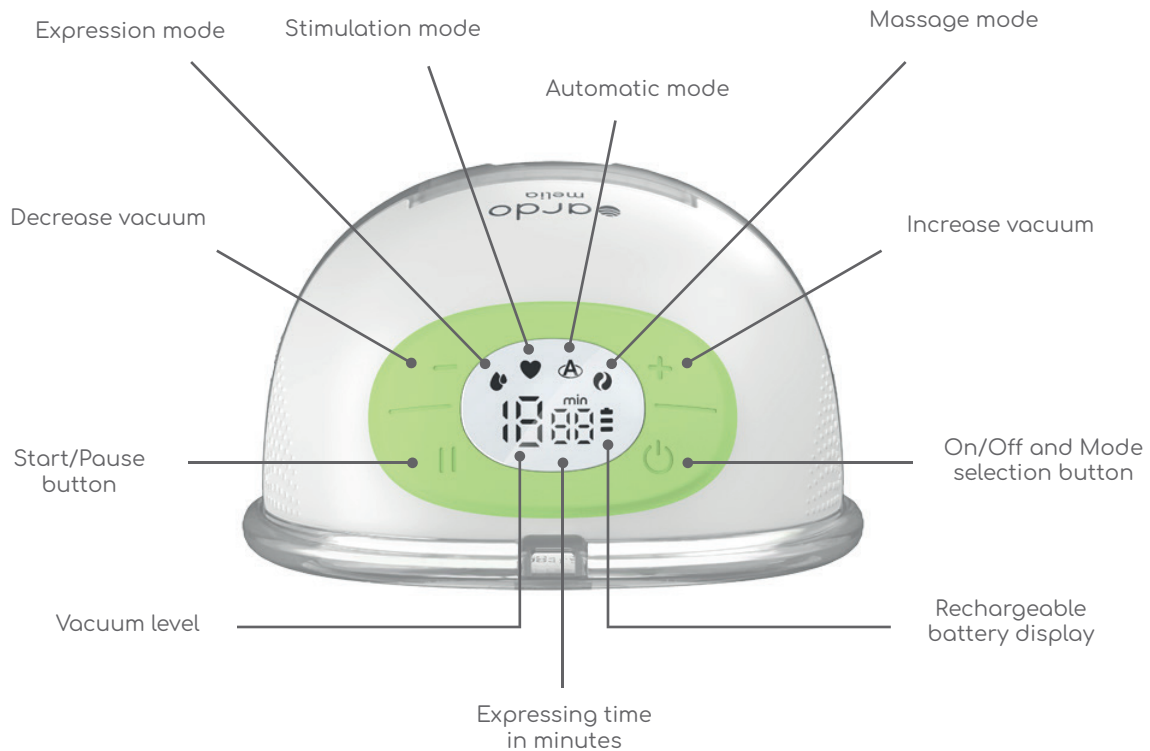

Intelligent LED display

- Clear information about the selected mode and vacuum strength
- Enables easy adjustment and monitoring during expression
- Can also be operated in the dark
- With pause button for short interruptions
- Automatic safety switch-off after 20 minutes

Easy handling

- The assembled breastpump can be placed upright for easy handling, as can the milk collection tray when separated from the pump motor.
- Pump motor with knobbed structure for easy grip.
- Quick and easy assembly and cleaning of the few individual parts.
- Except the pump motor: Clean in the dishwasher, sterilise in the microwave - practical and time-saving.