



KY-45SB
Basic type



KY-40QR
Emergency type



KY-45MB
Medicine type

Specifications

Product name		Medical cart		
Product code		KY-45SB	KY-40QR	KY-45MB
Type		Basic	Emergency	Medicine
Overall dimensions	Width	790 mm	900 mm	790 mm
	Length	518 mm	538 mm	543 mm
	Height	950 mm	1300 – 1830 mm (Including I.V. hanger rod)	1490 mm
Small-sized drawer		570 mm (W) x 395 mm (L) x 60 mm (H)		
Middle-sized drawer		570 mm (W) x 395 mm (L) x 135 mm (H)		
Large-sized drawer		570 mm (W) x 395 mm (L) x 215 mm (H)		
Casters		125mm with the diagonal locking system and one of the casters is the anti-static caster		
Material		Body: Powder-coated steel Cover: ABS plastic		
Coating		Powder coating		

Medical Cart

KY- 40 Series

- Supporting the management for medicine & emergency tools -



Emergency
use



Medicine
management

Basic type KY-45SB

Meets basic management needs for use in medical facilities

Only one drawer system with dividers

For safety reasons, the drawers are in the unlocked position, more than 2 drawers cannot be opened at the same time. Moreover, specially designed dividers are provided with 3 different sizes of drawers for orderly management of medicine and diagnostic materials.


Only one drawer system



Example of use



Additional dividers

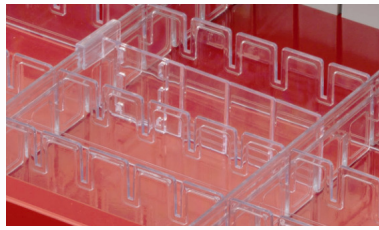

Always only one drawer is available for safety.

Handle



The cart comes with an ergonomically designed handle for easy transportation.

Push-style sub-table



A smart sub-table is available as needed depending on the type of work to be done in a hospital ward. The edges of the table are raised to prevent items from falling off. The table comes with small sections into which tools can be organized effectively.

Caster

Exclusive 125-mm casters have a diagonal-locking system to prevent the cart from being moved inadvertently. In addition, one of the casters is anti-static to ensure reliable use in a hospital ward.



Fully locking key mechanism

All drawers are fully locked with only one keyhole being locked. This function dramatically enhances work efficiency and prevents a drawer from being left unlocked.



Medicine type KY-45MB

Comes with added features for the “Management of medicine”



Transparent medicine cases with name plate

The specially designed transparent cases are equipped for the efficient management of medicines. The contents can be taken out immediately, and its transparent feature makes users easy to confirm the contents of cases.

Emergency type KY-40QR

Comes with added features for the “Use for emergency”



Ready for use state (Management stage)



Available for use state (Actual use stage)




Oxygen cylinder holder

The oxygen cylinder holder is provided as a standard accessory for the emergency cart. Note: A 500-liter-capacity oxygen cylinder is available.



Standard accessories for the Emergency type -IV hanger rod

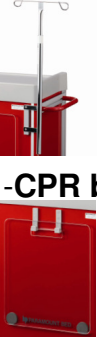

-CPR board with storage